PACKAGING


MEDIMASK N95

Protection Respirator

BOX - 30 PIECES

INDIVIDUAL BAG - 1 PIECE

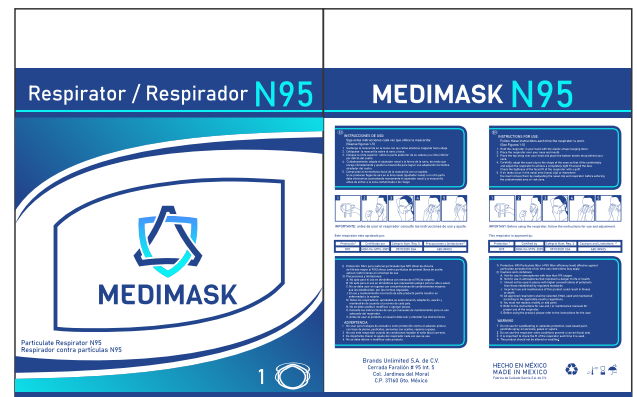BOX

MEDIMASK N95

Protection Respirator

www.medimask.mx

BAG

MEDIMASK N95

Protection Respirator

Respirator / Respirador N95

MEDIMASK

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INSTRUCTIONS FOR USE:
Follow these instructions each time the respirator is worn.
(See Figures 1-5)

1. Hold the respirator in your hand with the elastic straps hanging down.
2. Place the respirator over your nose and mouth.
3. Place the top strap over your head and place the bottom elastic strap behind your neck.
4. Carefully, adjust the nasal clip to the shape of the nose so that it fits comfortably and adjust the respirator to achieve a completely tight fit around the face.
5. Check the tightness of the facial fit of the respirator with a puff.

If air leaks occur in the nasal area (nasal clip) or elsewhere, you must remove them by readjusting the nasal clip and respirator before entering the contaminated area or risk zone.

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IMPORTANT: Before using the respirator, follow the instructions for use and adjustment.

This respirator is approved by:

Protection ¹	Certified by	Categoria Num. Reg. S.	Cautions and Limitations ²
N95	NOM-16-STPS-2009	093IC2020 SSA	ABC, JMMMS

1) Protection: N95 Particulate filter (95% filter efficiency level) effective against particulate aerosols free of oil. Use restrictions may apply.

2) Cautions and Limitations:

- A. Not for use in atmospheres with less than 10% oxygen.
- B. Not for use in atmospheres that represent a danger to life or health.
- C. Should not be used in places with higher concentrations of pollutants than those established by regulated standards.
- D. Incorrect use and maintenance of this product could result in illness or death.
- E. All approved respirators shall be selected, fitted, used and maintained according to the applicable country regulations.
- F. No must not clean, modify or add parts.
- G. Refer to the instructions for use and / or maintenance manuals for correct use of the respirator.
- H. Before using the product, please refer to the instructions for the user.

WARNING:

1. Do not use for sandblasting or asbestos protection, lead-based paint, pesticide sprays, oil aerosols, vapors or gases.
2. Do not use this respirator when conditions prevent a correct facial seal.
3. It is important to check the fit of the respirator each time it is used.
4. This product should not be altered or modified.

Brands Unlimited S.A. de C.V.
Cerrada Farallón # 95 Int. 5
Col. Jardines del Moral
C.P. 37160 Gto. México

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Fabrica de Calzados García S.A. de C.V.